


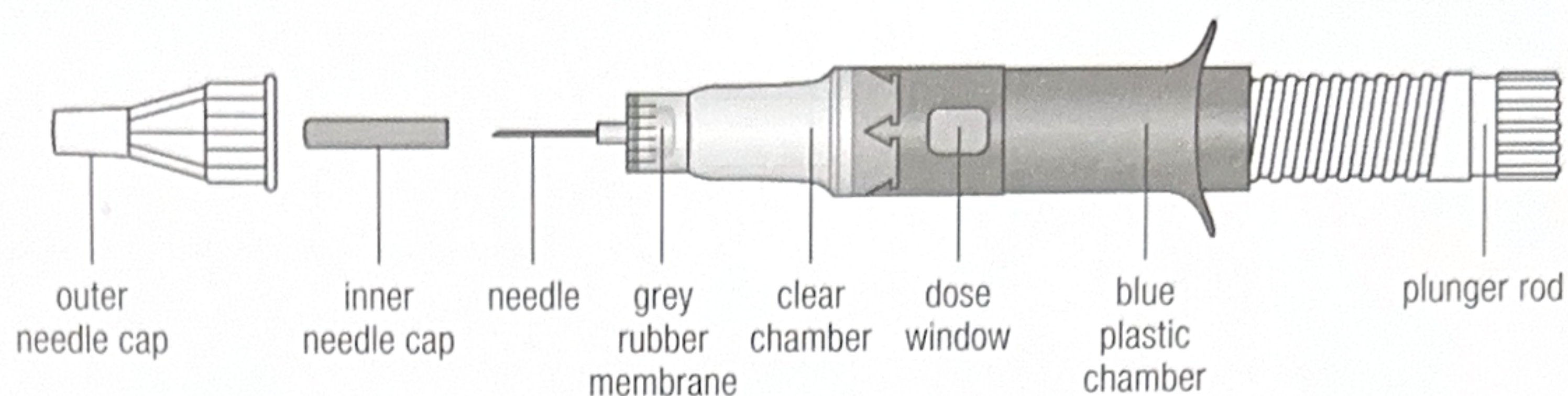


Before you use CAVERJECT[®] Impulse. My prescribed dose is _____ micrograms.

- Wash and dry your hands before using CAVERJECT[®] Impulse. If your partner wishes to be involved in preparing or giving the injection, they should also wash and dry their hands.
- CAVERJECT[®] Impulse is a **single use** device which can administer only one dose.
- It is important that you follow the instructions carefully to prepare and inject CAVERJECT[®] Impulse.
- It is important that you also read the Consumer Medicine Information leaflet available from your pharmacist and follow your doctor's advice.


You can prepare and store your CAVERJECT[®] Impulse for up to 24 hours before using it. Prepared Caverject Impulse must be stored in the refrigerator at 2-8°C. Do Not Freeze.

How to use CAVERJECT[®] Impulse

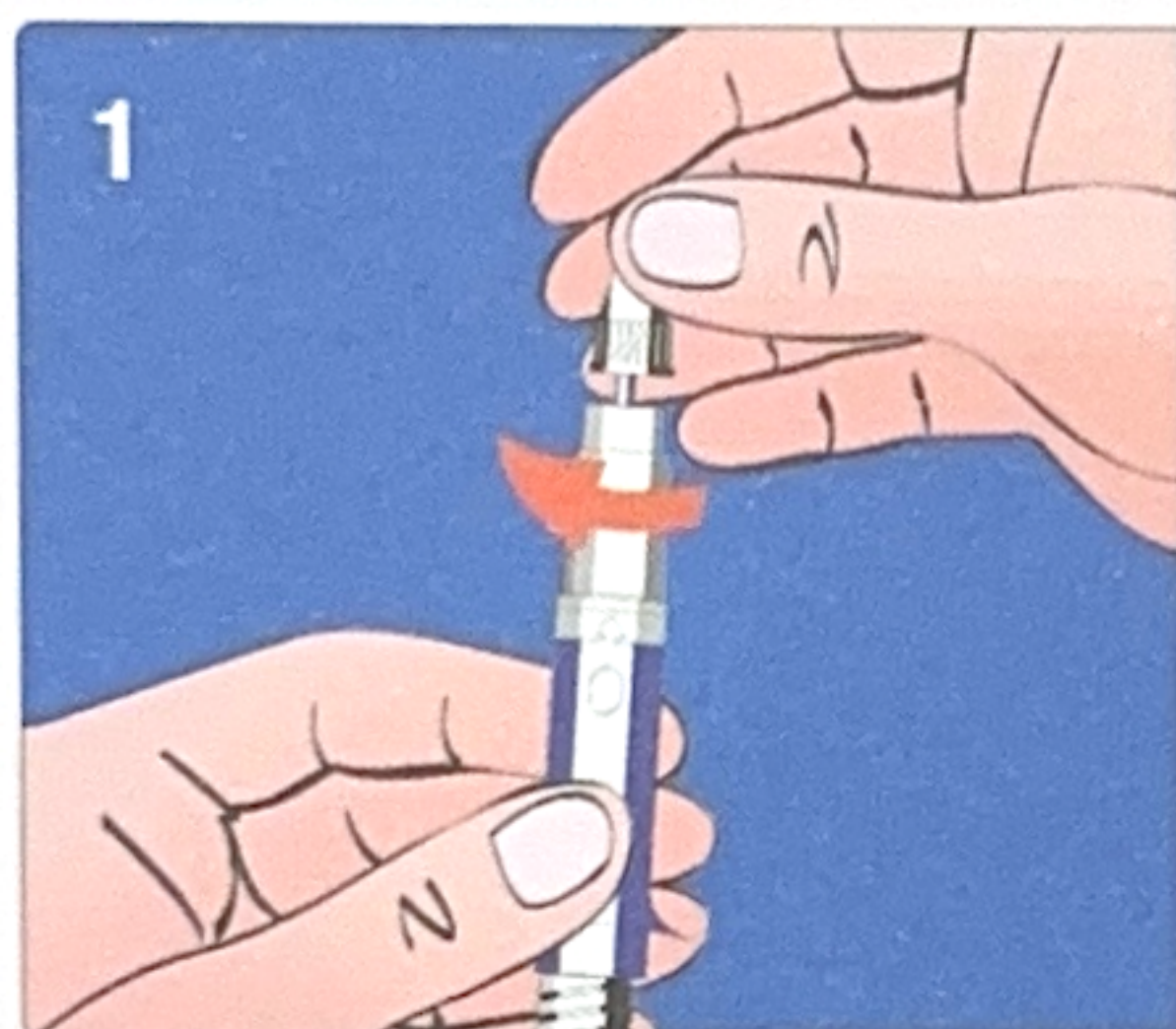


During this step you will be connecting the needle to the device.

- Remove device and capped needle from the box.
- Clean the grey rubber membrane at the tip of the device with one of the alcohol swabs provided.
- Completely peel the foil from the needle cap.
- Hold the syringe with the needle-end pointing upwards.
- Attach the needle onto the tip of the device by pushing and turning until it is firmly in place.

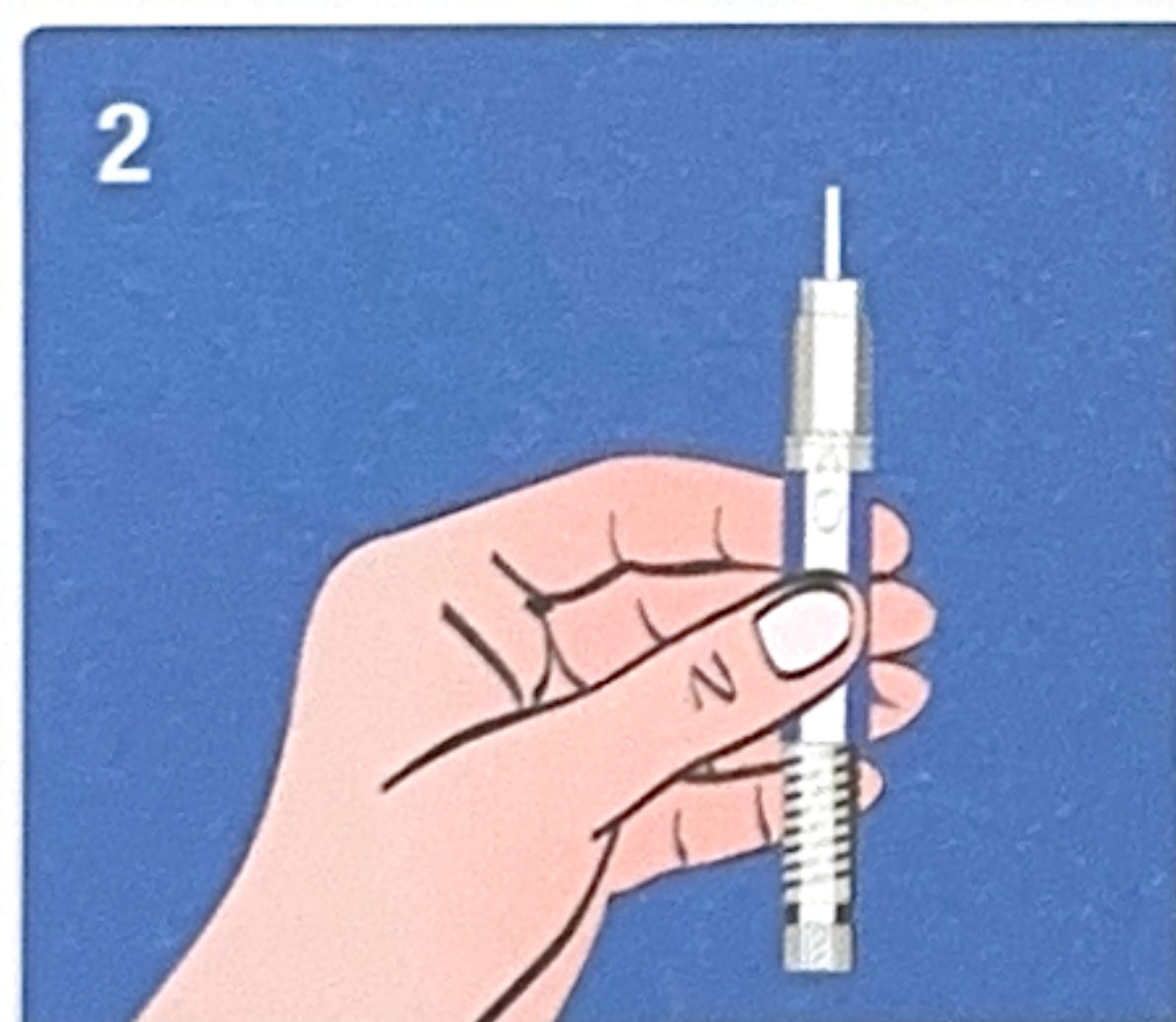


During this step you will be removing the outer needle cap.

- Remove the outer cap from the needle.
- Hold the syringe with the needle-end pointing upwards.
- The plunger rod will be in the extended position.

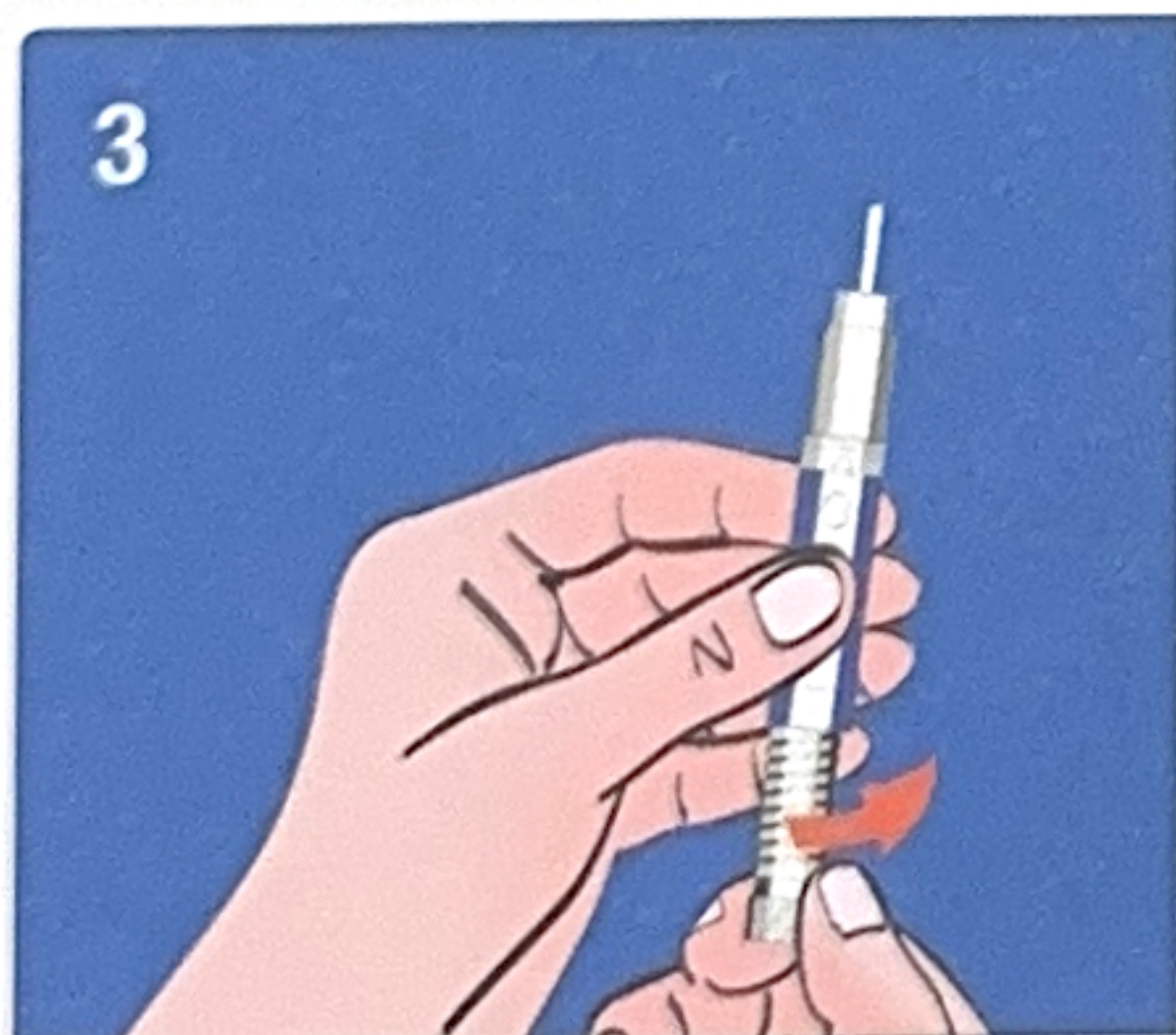

During this step you will be mixing the powder and the liquid in the device.


- Hold the device with the needle-end pointing upwards, gripping the blue chamber with one hand.
- Turn the plunger rod as shown left until it will not turn any further (fully wound).
- The reconstituted solution will now appear in the clear chamber.
- The solution may appear cloudy at first. Gently shake to make sure the solution is evenly mixed and the powder has dissolved. Do not use if the solution remains cloudy or contains particles.

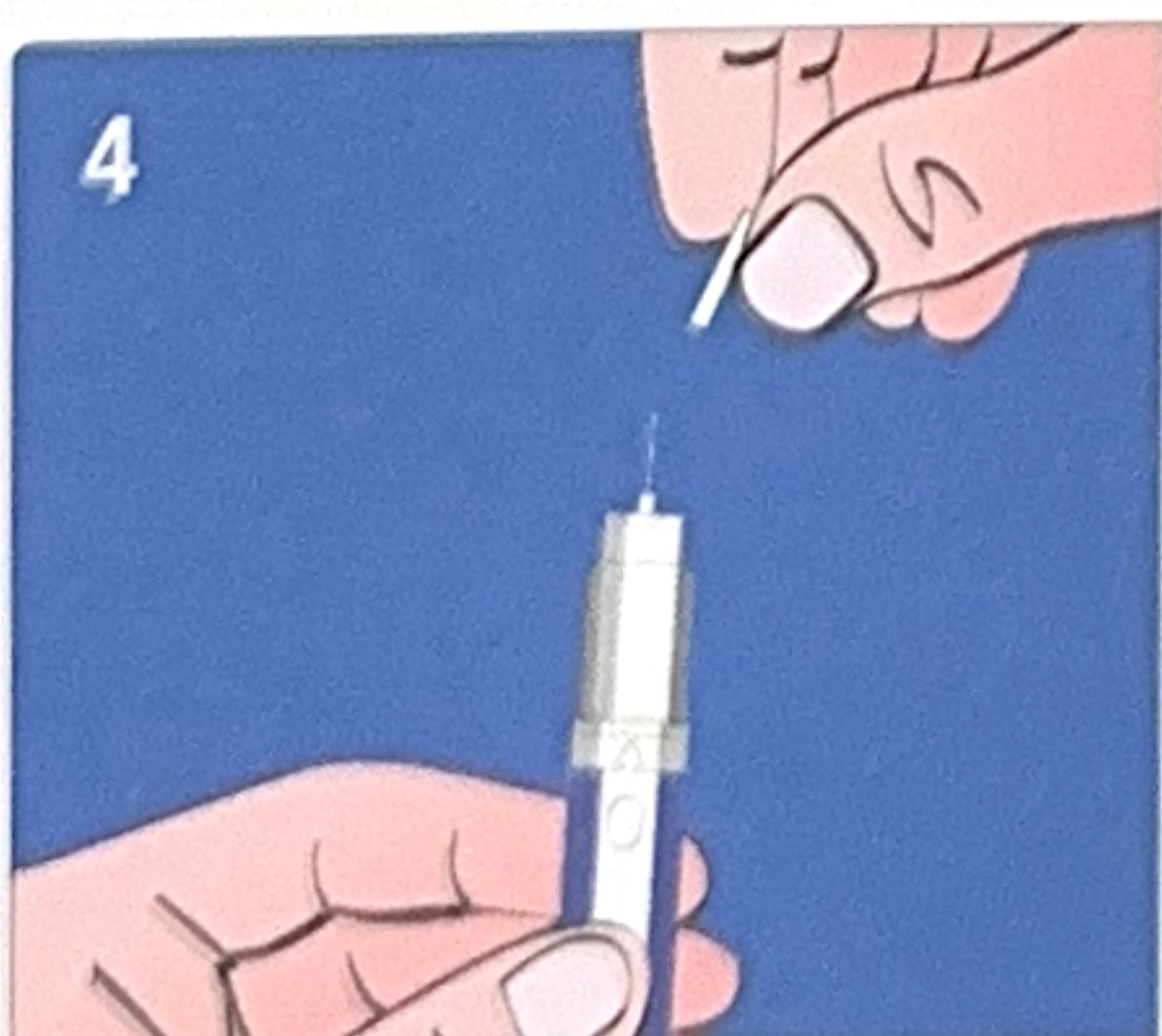

During this step you will be removing the inner needle cap.

- Hold the device with the needle pointing upwards.
- Carefully remove the inner cap from the needle.
- Do not touch the tip of the exposed needle.
- The needle is very fine. Take care not to bump it.
- **Do NOT use if the needle is bent.**

Continued over

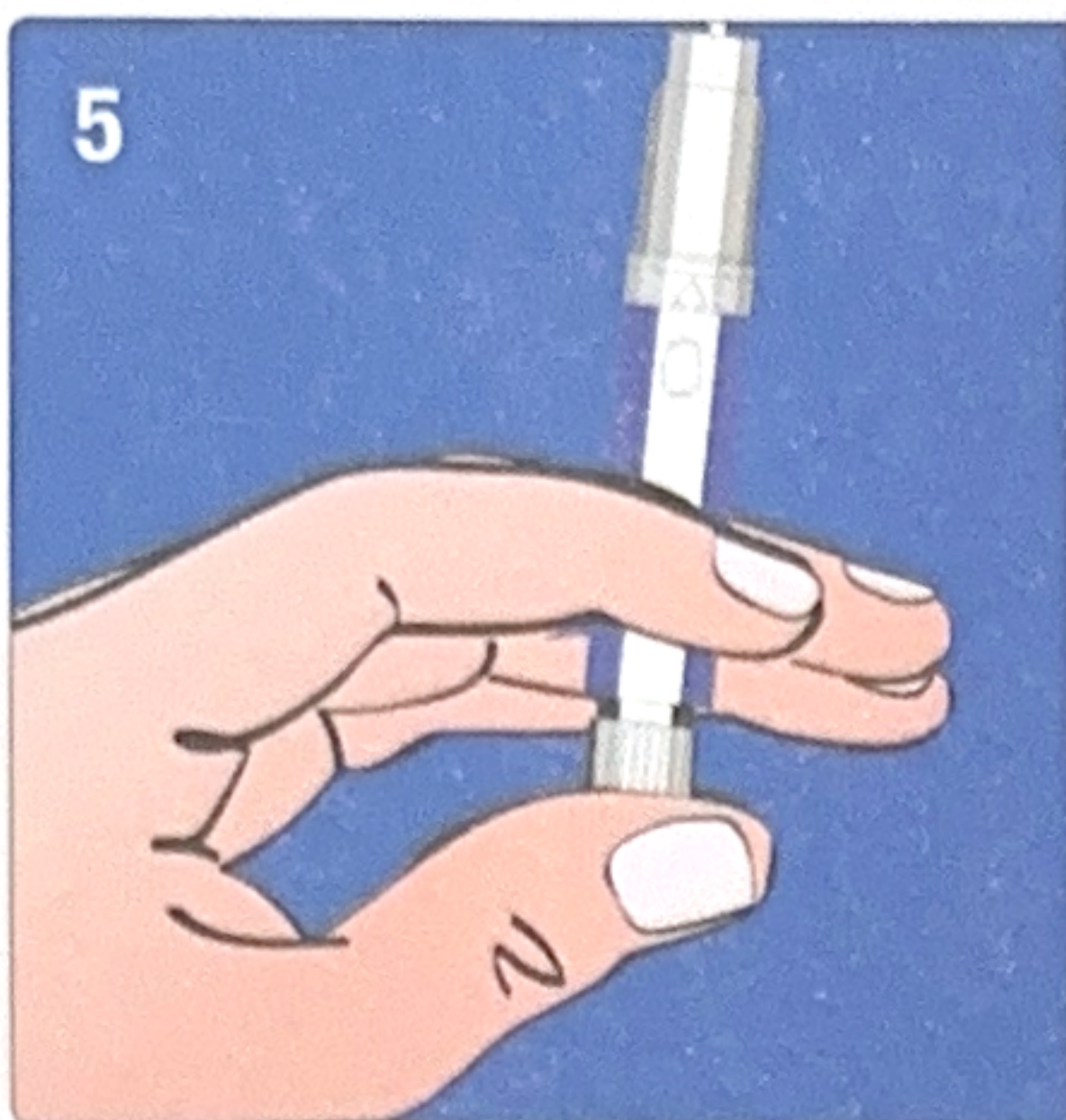

5

During this step you will be removing air from the device. A few drops will appear at the needle tip. This is expected.

- Remove any large bubbles by gently tapping the clear chamber with your finger.
- Ensure the plunger rod is still fully wound.
- Hold the device in **one hand** as shown, with needle-end pointing upwards.
- Remove the air from the clear chamber by pressing the plunger rod with the thumb as far as it will go.


If you need to use both hands to complete this step, the plunger rod must be pushed directly upwards, **do not twist**.



6 Your prescribed dose #

During this step you will be dialling your dose.

- Hold the device in one hand.
- Grip the blue chamber with the dose window facing you (needle tip to the left).
- Turn the plunger rod slowly.
- Dial your prescribed dose.
- If you make a mistake, keep turning the plunger rod until the correct dose reappears.
- Set down the device on a level surface ensuring needle tip does not make contact with anything.

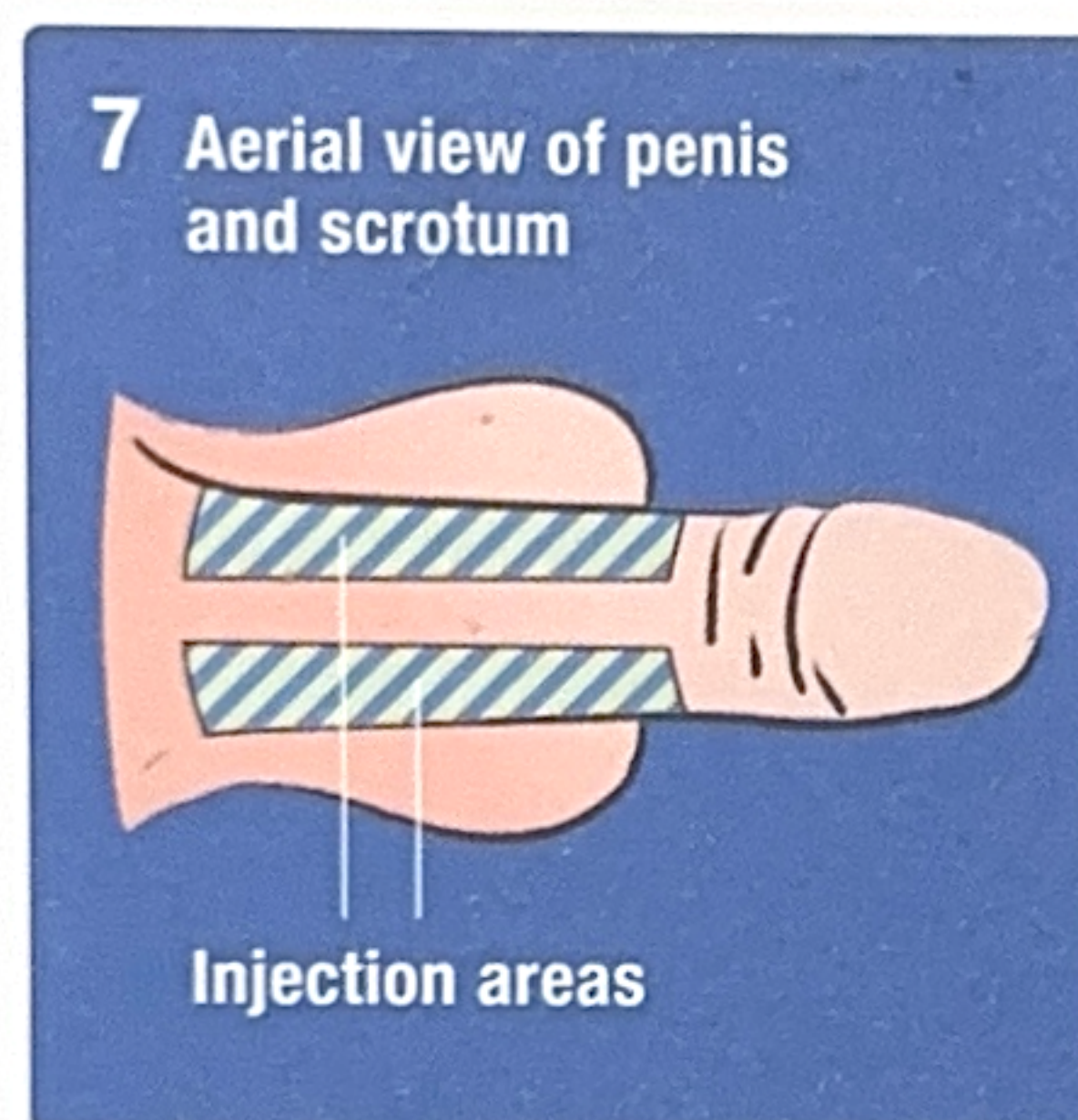


7 Aerial view of penis and scrotum

In this step you will be preparing the injection site on your penis.

CAVERJECT® Impulse is injected once into the side of your penis on each occasion you use it. It is recommended to alternate the site of injection on each occasion. See diagram left.

- Sit in an upright or slightly reclined position.
- If you are uncircumcised, pull back the foreskin.
- Hold the head of the penis straight ahead.
- Clean the site selected with an alcohol swab.

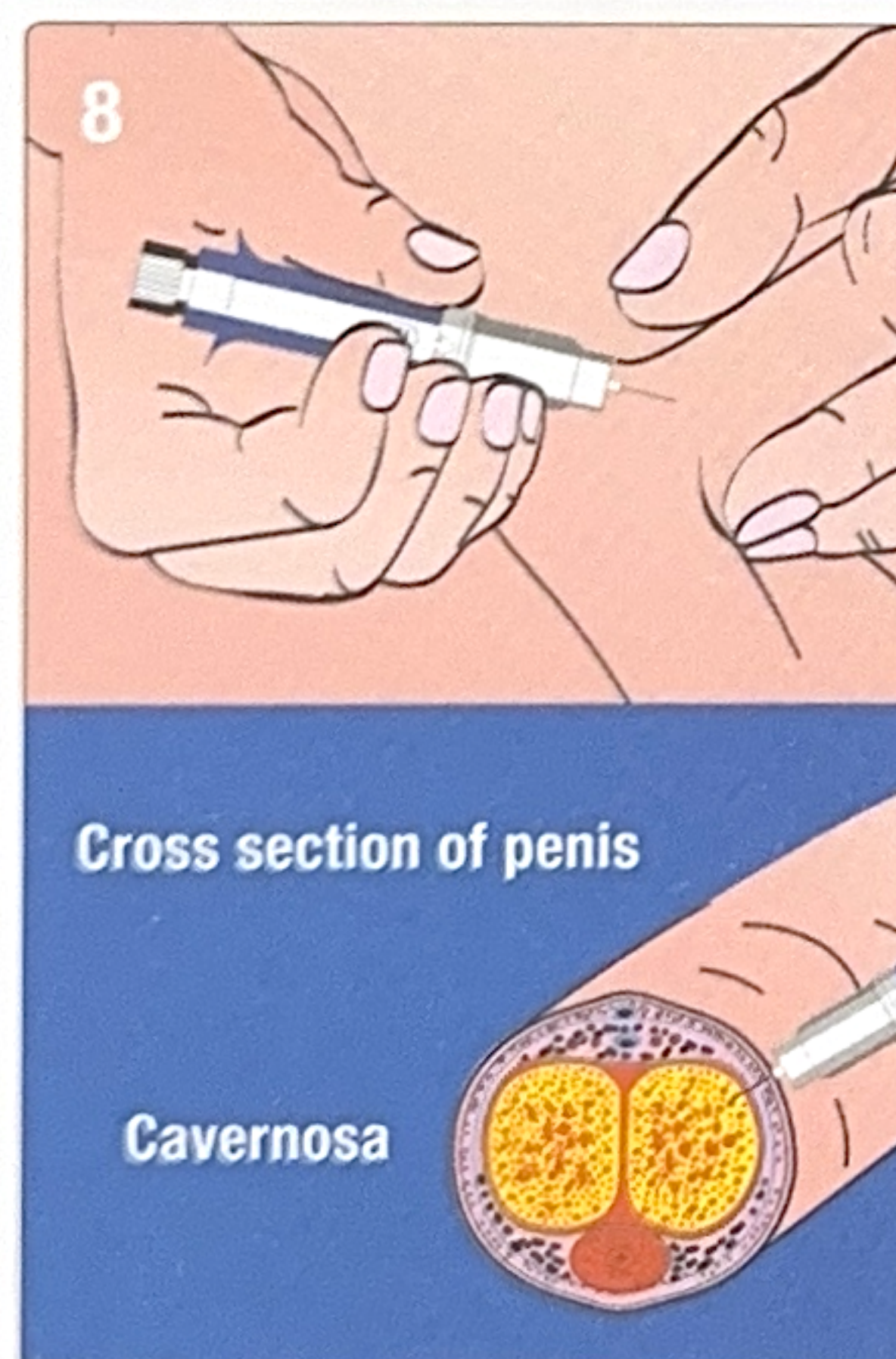

8

In this step you will be inserting the needle and injecting the dose into your penis. The injection needs to go into the cavernosa of the penis as shown in diagram left. Ensure you inject into the side of your penis and NOT into the top or underside. Avoid any visible veins.

- Gently rest the needle tip on the side of your penis where you have cleaned.
- Ensure the dose window is facing you and your dose is still centred.
- The needle should be inserted at an angle of 90°. With a steady and continuous motion insert the entire length of the needle.
- Press the plunger rod with your thumb. Injecting the dose requires you to apply some pressure to the plunger rod. This is normal.
- Withdraw the needle and apply pressure to the injection site with the alcohol swab for about 5 minutes or until any bleeding stops.
- Massage the penis to help CAVERJECT® Impulse spread through.
- Depending on your dose some solution may remain in the device. The device can be used only once, after which it should be disposed of safely.

For further information on the correct use of CAVERJECT® Impulse, please call 1800 675 229 or visit www.cject.com.au

REMEMBER: If you have any doubts or questions, ask your doctor for advice.

